

How to perform the Ling sound check six sounds + silence

- If the child wears two hearing devices, put both of them on.
- Make sure both devices are turned on.
- Remove the coil or hearing aid at the right ear so the child is using only the left device.
- Place yourself next to or behind the child so she can hear you, but can't see your face.
- Say the following sounds one at a time and teach the child to repeat each sound after you say it: "ah," "oo," "ee," "s," "sh," "m." Include a moment when you say nothing (silence).

ahhh



ooo




eee



shhh



sss

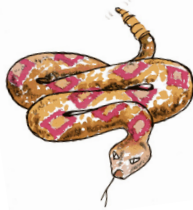

mmm



silence



- Each time you perform this check, say the Ling sounds in a different order. Otherwise, the child might memorize the sounds one way and just repeat them back, without actually listening for a prompt.
- Use unpredictable intervals of silence between each sound so the child has to respond only after she hears you say something.
- Replace the coil or device at the right ear and remove the coil or device from the left ear so the child is using only the right device.
- Repeat the Ling check with that device.
- Replace the coil or hearing aid so the child is now using both devices.
- Note any difficulty the child has repeating these sounds back to you with each device.